Professional
SHIATSU-AKUPRESSUR

Art. 0345

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Dear Customer,

Thank you for purchasing one of our electronic operated Shiatsu acupressure devices. Quality and design are major factors in a world that is controlled by mass production. Our brand Peter Bausch stands for high quality. We associate our name and the composition of our products with excellence and offer products which provide pleasure and enjoyment. The satisfaction of our customers is most important to us!

SHIATSU ACUPRESSURE SET ITEM 0345

This set with 45 W power engine, includes a mains adapter (adjustable to 5 speeds), 4 different attachments for diverse needs (a – pointed flat 22 mm, b – half pointed 22 mm, c – pointed 20 mm, d – rounded flat 25 mm) and a long spiral cord for free moving space.

1. Location hole
2. Speed regulator
3. Adapter
4. Small plug
5. On / Off
6. Pilot lamp
7. Spiral cabel

 Please read all available documentation (instruction of use and safety instructions) carefully before using the product! We assume neither responsibility nor liability for accidents, damages and malfunctions resulting from inobservance of these guidelines. The instructions have to be stored in a safe place and have to accompany the product at all times; Especially, if the product gets passed on to a third party!

 Never use the device in the bathtub, shower or if sink is filled with water! Never clean the device with water! Always use a damp cloth for cleaning!

- Unplug after usage!
- Children should never use this device without direct supervision of an adult!
- Always check the device before using it! If there are any visible damages, DO NOT USE IT!
- Attention! Do not open or damage the device by forcefully impacts. Injury and electrical shock hazard!

SAFETY INSTRUCTIONS

- Only plug in appropriate isolated ground receptacles (220-240V)!
- Only usable with the enclosed mains adapter (Item 0239)!
- Do not bend, crush or squeeze the connecting cable! Risk of electrical shock due to damaged cable insulation!
- Connecting cable should never come in contact with hot surfaces! Risk of electrical shock due to damaged cable insulation!

STARTUP

1. Press attachment tightly into device (1). Insertion will be confirmed by a "click"-sound.
2. Connect the device to the mains adapter (4). Pay attention to the proper polarity (+/-). If everything is linked together correctly, the control lamp (6) will light up. (Note: The revolution speed will be adjusted on the mains adapter surface.)
3. Adjust revolution speed: Slow speed: level 1 = 3 V / level 2 = 4,5 V
Average speed: level 3 = 6 V / level 4 = 7,5 V
High speed: level 5 = 9 V
4. Turn device on (5).
5. To protect the mains adapter: Never use the device longer than 15 minutes.

GENERAL INFORMATION

Massages are the natural way of gaining more vitality and wellness in your life. Your blood and lymphatic circulation will improve considerably through massages. Slow and extensive massages affect muscle relaxation. Whereas punctual and deep massages cause an enhancement of body tissue tension due to the stimulation of arterial blood and lymph flow. The exchange of body fluids as well as the nutrition of tissue will be supported. A quicker regeneration of injured muscle tissue is ensured by oxygen in the blood and the stronger blood circulation boots the metabolism effectively. Furthermore, muscular fatigue will be reduced and the musculature will be relaxed.

APPLICATION

Besides the stimulation of individual organs, the device can be used for whole body massages and reflex zone therapy. A continuous appliance of the BAUSCH Shiatsu acupressure device is necessary for the desired success. Pay attention to your body's reaction. Avoid any exaggerated implementation. The following illustrated charts will give you an overview about the attachments and their field of application as well as the recommended rotation speed.

Consult a doctor, if you have any doubts about your health condition before using the device!

CLEANING & MAINTENANCE

Always use a damp cloth for cleaning the device!
Never rinse it under water!

SPECIAL INSTRUCTIONS

The high-quality BAUSCH Shiatsu acupressure device is equipped with a complex technology which enables it to convert the rotary motion from the engine in a low friction "walk movement". Because of the reason given a mechanical noise development is inevitable due to the function of multiple long-lasting metal bearings.

An excess of oil residue under the attachments might be possible after the first couple of applications. This does not endanger your health at all! Please clean it with a soft cloth.

REPAIR INFORMATION

Electric devices are only to be fixed by qualified personnel. Improperly repairing can result in severe consequential damage.

CAUTION!

Never use the device if you suffer from acute inflammation (e.g. nephritis, phlebitis or arthritis); acute tendon, ligament, nerve or muscle injuries; prolapsed intervertebral disc; bone diseases like osteoporosis; diseases of the vascular system like thrombosis; fever; cancer; heart diseases; frostbite or burns.

Precaution if you suffer from skin cancer (including inflammations under the skin) and varicose veins.
If you are pregnant or suffer from diabetes we advise consulting a doctor.

DISPOSAL

 This product as well as any other electronic appliances, which have reached the end of their lifetime, must be collected and disposed of separately from household waste material. You can get information on where to dispose your old devices free of charge from your community or municipal administration.

IF YOU HAVE ANY QUESTIONS, FEEL FREE TO CONTACT US:
info@peter-bausch.de

Stufe	1	2	3	4	5
Volt	3	4,5	7	7,5	9

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I. STIMULATION OF INDIVIDUAL ORGANS/AREAS


Level	1	2	3	4	5
Volt	3	4,5	7	7,5	9

Intention of the Treatment	Area of Application		a	b	c	d	Comments
Circulation of the head, arms and hands. Also migraine prophylaxis.						5	Do not use during migraines
Improvement of blood pressure.						4	Smooth out heart area at the back
Easing of heart nerves and the heart itself.		Also foot reflexology				4	Slight pressure stimulates the nerves. Heavy pressure calms the nerves
Mood-lifting (if you suffer from anxiety).						4	Heavy vibrations on costal arch
Relief of limp pain.		Also foot reflexology			3		Use pressure until pain ceases
Relief of upset stomach.		Also foot reflexology		3			Slightly smooth out the stomach area. Do not use if you suffer from inflammatory disease
Support of the digestive organs.		Also foot reflexology		3			Slightly smooth out the stomach area. Do not use if you suffer from inflammatory disease

Intention of the Treatment	Area of Application		a	b	c	d	Comments
Support of the intervertebral discs.						5	Smooth out gently, gentle vibrations. Never treat the spine directly. Only if there are no medical concerns!
Reduction of blood pressure.				4			Smooth out the whole back
Tightening up the abdominal muscles.						4	Massage heavily until redness is visible. Only if there are no medical concerns!
Relief of insomnia.						5	Massage calves heavily. Smooth out the back
Reduction of cellulite.						4	Heavy massages preserve an elastic tissue
Circulation of the legs, reduction of swellings.						5	Massage heavily until redness is visible
Relief of sciatica troubles.		Also foot reflexology		4			Never massage if acute! Do not use if you suffer from kidney diseases!
Tightening up the gluteal muscles.						4	Massage heavily until redness is visible

Attention: If you have a mains adapter that adjustable to 6 speeds: DO NOT USE LEVEL 6 (12V)! The impact on your skin will be too heavy!

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
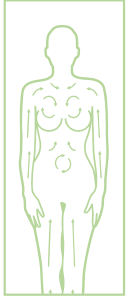
3.1 FUSSREFLEXZONENTHERAPIE

3.2 HANDREFLEXZONENTHERAPIE

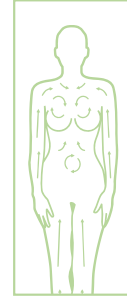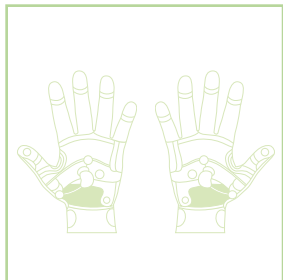The entire organism is positively influenced by the effectiveness of the whole body massage. Basically: Skin and muscle tissue is supplied with blood perfectly.

Whole body massages always concentrate on massaging in direction to the heart. Start from the right foot and work your way up the right leg. Repeat this step for the left foot and leg. Now continue the massage: first do your right hand and arm, then switch to the left side and massage your hand and arm. Proceed with your upper body part.

The massage concludes in the area of the heart. Do not massage dorsal process; however, massaging the lateral muscle areas can be very affective!

3. REFLEX ZONE THERAPY

Foot and hand reflexology deals with the assumption that different reflex zones on the feet and hands can stimulate various body parts and organs.


By massaging those zones it is possible to influence the appropriate inner organs and body parts. This massage is suitable for regular health care and the treatment of smaller discomforts.

A reflex zone therapy comprises of putting pressure on a specific point on hands or feet. The BAUSCH Shiatsu acupressure device and its various attachments enhance this procedure. The massaging head compensates the finger and functions as the contact point for the pressure massage. All of the reflex zones can be massaged ideally due to the various types of attachments. There is no risk of injuries caused by fingernails. Furthermore, a considerably more burden-free massage is ensured.


Massage the chosen point with one of the attachments. Either use slight or heavy pressure on this point to self-regulate the massage's effect.

Moreover, you can regulate the speed of the device. By choosing a higher speed level the mechanical stimulus will be raised also. Remember: The massage should never exceed your level of pain. After multiple appliances, your level of pain will have been decreased which is an indicator for a successful therapy.


3.1 FOOT REFLEXOLOGY

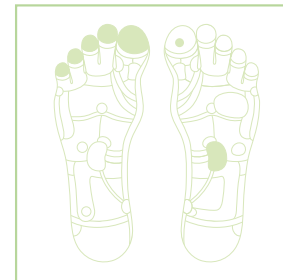
3.2 HAND REFLEXOLOGY

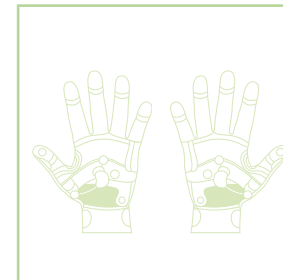Stufe/Level	1	2	3	4	5
Volt	3	4,5	7	7,5	9

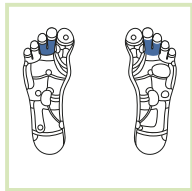
Augen
Eyes
c / 4

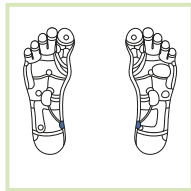
Blase
Bladder
c / 5


Harnleiter
Ureter
c / 3-4

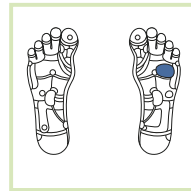
Herz
Heart
b / 3-4

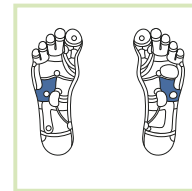
Leber / Gallenblase
Liver / gall bladder
d / 5

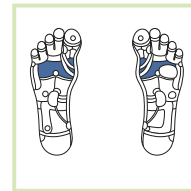
Lungen / Bronchien
Lungs / bronchial tubes
d / 5

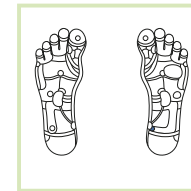
Mastdarm
Rectum
c / 4

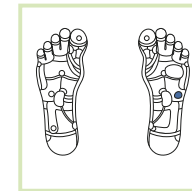
Milz
Spleen
b / 4

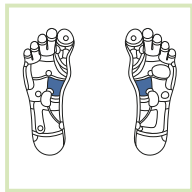
Magen
Stomach
d / 4

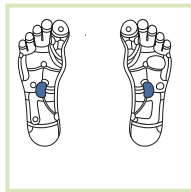Nieren
Kidneys
a / 4-5

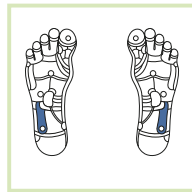
Dickdarm, absteigender
Large intestine, descending
b / 5

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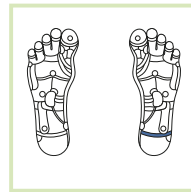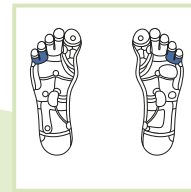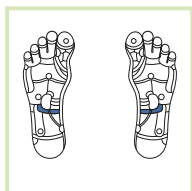
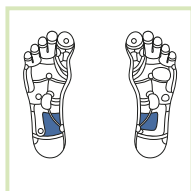
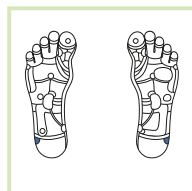
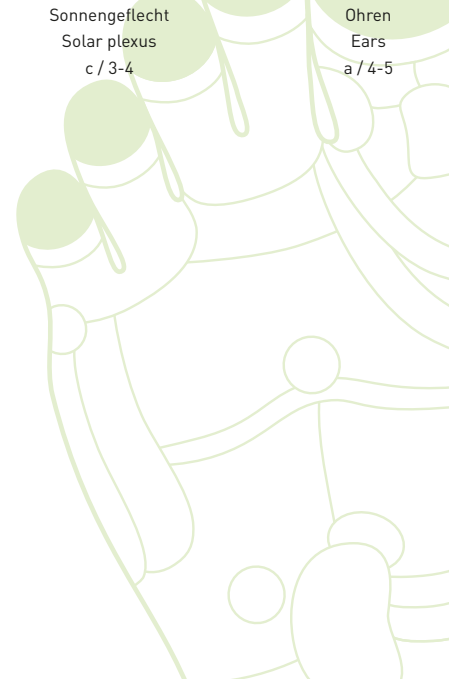
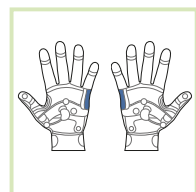
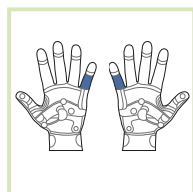
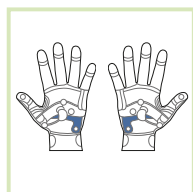
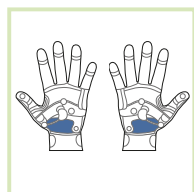
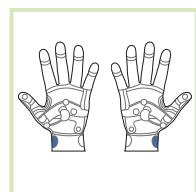
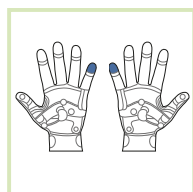
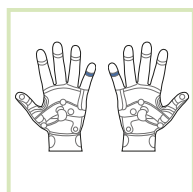
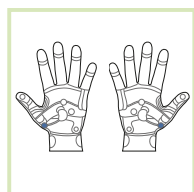
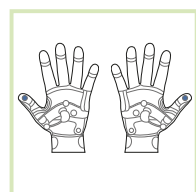
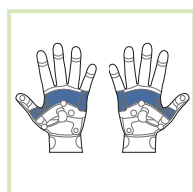
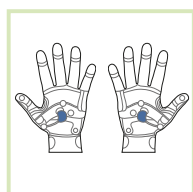
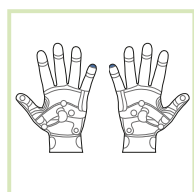
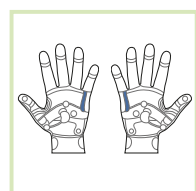
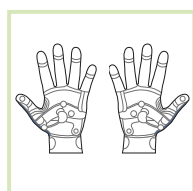
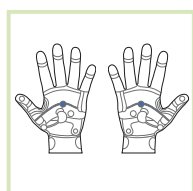
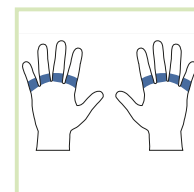
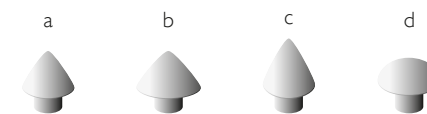
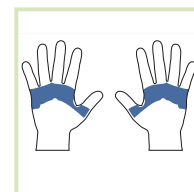
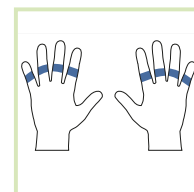
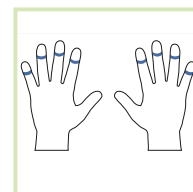
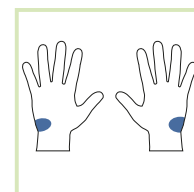
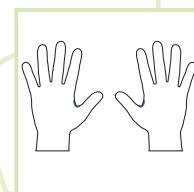
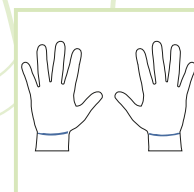
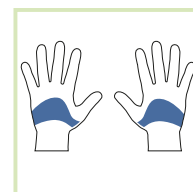
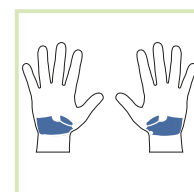
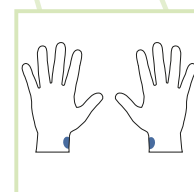
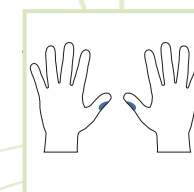
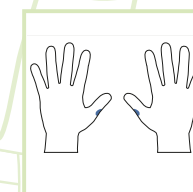
■ **CE-MARKING**

The CE-marking testifies that the relevant, essential requirements of the appliance comply with the pertaining European norms.